

The CONTACT Wax Dispenser does not only impress with the elegant design, it offers equally great benefits for everyday work. The new Contact Wax Dispenser combines aesthetics and functionality and allows the waxes to adhere in place individually. The slight inclination of the plate allows clear and ergonomic work. A small useful detail is the practical wax scraper for cleaning the modeling tips.


excellence
line

CONTACT



CONTACT Wax dispenser

CONTACT = Contact Basis =
Base plate with wax scraper and dust cover

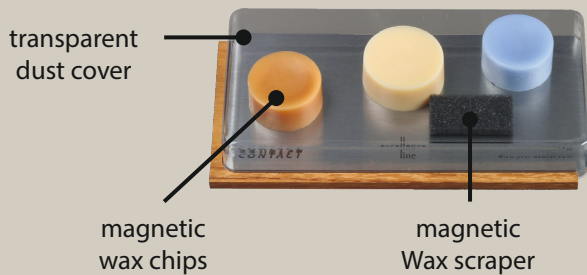

optional accessories:



CONTACT-E



CONTACT-E hot Wax dispenser

- **Contact base** (see CONTACT)
- **the left half of the metal plate is heated**
- fixed temperature of approx. 35 - 40°C
- can be switched on and off

CONTACT-E Wax Up Station


- **Contact base** (see CONTACT)
- **the left half of the metal plate is heated**
- fixed temperature of approx. 35 - 40°C
- can be switched on and off
- **incl. 1 handpiece with modeling tip PKT2**
- infinitely variable temperature control up to 190°C
- incl. handpiece clamp for safe deposit / clamping of the handpiece after use

CONTACT-E duo Wax Up Station

- **Contact base** (see CONTACT)
- **the left half of the metal plate is heated**
- fixed temperature of approx. 35 - 40°C
- can be switched on and off
- **incl. 2 handpieces with modeling tip PKT1 and PKT2**
- separate stepless temperature control up to 190°C
- incl. 2 handpiece clamps for safe depositing/clamping of the handpiece after use

The unique mobile phone holder also adheres magnetically to the CONTACT. In the future, helpful tutorials "wax up" can be accessed directly at the workplace via Internet. Acquired knowledge can be implemented immediately "step by step". To protect against contamination, the transparent dust cover can simply put over after the work.

„CONTACT“ Wax dispenser / „CONTACT-E“ Wax Up Station overview / refill packs

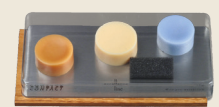Contact Basis = Base plate with wax scraper and dust cover			
	780-0100 CONTACT freestyle Wax dispenser - Contact Basis		739-6320 Effect - wax Chip CONTACT blue transparent / orange transparent / incisal transparent (ashfree) each 20g
	416-0000 CONTACT - E hot Wax dispenser - Contact base - left side tempered		738-6325 Sculpturing wax Chip CONTACT sahara (ashfree) by J. Peters, 3x25g
	416-1000 CONTACT - E Wax Up Station - Contact base - left side tempered - with 1 handpiece + 1 modeling tip (PKT2)		738-6335 Sculpturing wax Chip CONTACT ivory (ashfree) by J. Peters, 3x25g
	416-2000 CONTACT - E duo Wax Up Station - Contact base - left side tempered - with 2 handpieces + 2 modeling tips (PKT 1, PKT2)		735-6325 Sculpturing wax Chip CONTACT beige (opak) 3x25g
      			734-6325 Sculpturing wax Chip CONTACT grey (opak) 3x25g
Art.-No. 411-.... 0100 0200 0300 0400 0500 0600			737-6325 Sculpturing wax Chip CONTACT emerald-green (ashfree) 3x25g
	780-0200 CONTACT traditional Wax dispenser Contact Basis, sculpturing wax chip blue (ashfree) 25g, sculpturing wax chip grey (opaque) 25g, Cervical wax chip red (ashfree) 20g		736-6325 Sculpturing wax Chip CONTACT blue (ashfree) 3x25g
	780-0300 CONTACT n. J. Peters Wax dispenser Contact Basis, sculpturing wax chip ivory (ashfree) 25g, sculpturing wax chip sahara (ashfree) 25g, Cervical wax chip wood (ashfree) 20g		733-6325 Sculpturing wax Chip CONTACT yellow (ashfree) 3x25g
	780-0400 CONTACT esthetik Wax dispenser Contact Basis, sculpturing wax chip blue transparent (ashfree) 20g, sculpturing wax chip orange transparent (ashfree) 20g, sculpturing wax chip sahara (ashfree) 25g		738-6320 Cervical wax chip CONTACT wood (ashfree) by J. Peters, 3x20g
	780-0500 CONTACT NAT Wax dispenser Contact Basis, sculpturing wax: blue, grey, neon-yellow, neon-green Underlining wax : caramel, Cervical wax: red 6 x 8g (ashfree) with magnetic tape		731-6320 Cervical wax chip CONTACT red (ashfree) 3x20g
			102-0002M IQ training set CONTACT by G. Seubert (ashfree), 6 x 8g with magnetic tape sculpturing wax: blue, grey, neon-yellow, neon-green Underlining wax : caramel, Cervical wax: red
			780-0010 phone holder CONTACT 1 pc
			780-0006 Wax scraper CONTACT 6 pcs